

NexiusZoom

HIGHLIGHTS

- ✓ Bino and Trinocular stereo zooms
- ✓ Zoom 1:6.7 ratio, 0.67x to 4.5x
- ✓ Ergonomically designed stands
- ✓ 3 W LED illuminators
- ✓ 10 years warranty



NexiusZoom with pillar stand
NZ.1902-P



NexiusZoom with rack and pinion stand
NZ.1903-S

Specifications

Eyepieces

Pair WF 10x/22 mm eyepieces

Head

Bino- or trinocular heads with 45° inclined tubes

Both eyepieces with ± 5 diopter adjustments

Interpupillary distance adjustable between 54 mm and 75 mm

Objective

Zoom 1:6.7 objective with 0.67x to 4.5x magnifications

Field of view from 33 mm to 4.9 mm

Working distance 110 mm

Auxiliary 0.5x / 0.75x / 1.5x and 2.0x lenses available

All optical parts are anti-fungus treated

Stands

Ergonomically designed pillar or rack and pinion stand with LED illumination and 2 object clamps

Ergonomically designed universal one arm or double arm stands without illumination

Alloy metal cast, hardened off-white painted

Focusing adjustment

Coarse adjustment

Illumination

The 3 W transmitted and incident LEDs illuminators can be used simultaneously

Both light intensities are adjusted separately

	Features					
Models	Binocular	Trinocular	Pillar stand	Rack & pinion stand	One arm stand	Double arm stand
NZ.1902-P	•		•			
NZ.1902-S	•			•		
NZ.1902-U	•				•	
NZ.1902-B	•					•
NZ.1903-P		•	•			
NZ.1903-S		•		•		
NZ.1903-U		•			•	
NZ.1903-B		•				•

Models



NZ.1902-P



NZ.1903-P



NZ.1902-S



NZ.1903-S




NexiusZoom with one arm stand
NZ.1902-U




NexiusZoom with double arm stand
NZ.1903-B

Accessories & Spare parts

NZ.5302	NexiusZoom binocular head with eyepieces
NZ.5303	NexiusZoom trinocular head with eyepieces
NZ.6010	Pair HWF 10x/22 mm eyepieces
NZ.6015	Pair HWF 15x/16 mm eyepieces
NZ.6020	Pair HWF 20x/12 mm eyepieces
NZ.6110	HWF 10x/22 mm eyepiece with micrometer
NZ.6099	Pair of eyeshades
NZ.8905	Auxiliary lens 0.5x, working distance 187 mm
NZ.8907	Auxiliary lens 0.75x, working distance 125 mm
NZ.8915	Auxiliary lens 1.5x, working distance 50 mm
NZ.8920	Auxiliary lens 2.0x, working distance 33 mm
NZ.8950	Protection glass for NexiusZoom head
NZ.9010	Ergo rack and pinion stand with transmitted and incident LED illuminators
NZ.9015	Ergo pillar stand with LED transmitted and incident illuminators
NZ.9020	Universal one arm stand without NexiusZoom head holder
NZ.9030	Universal double arm stand without NexiusZoom head holder
NZ.9090	NexiusZoom head holder for NZ.9020/9030
65.980	Articulated arm stand with table clamp without NexiusZoom head holder
68.981	Head 76 mm holder for black universal stand 65.980
NZ.9570	Pair object clamps for stage
NZ.9833	C-mount adapter with 0.33x lens for 1/3" cameras Only for NZ.1903 models
NZ.9850	C-mount adapter with 0.5x lens for 1/2" cameras Only for NZ.1903 models
LE.1973	Ring light illuminator with 144 LED's and segment controller
AE.1112	Object micrometer 50 mm divided in 500 parts on glass slide 76 x 26 mm

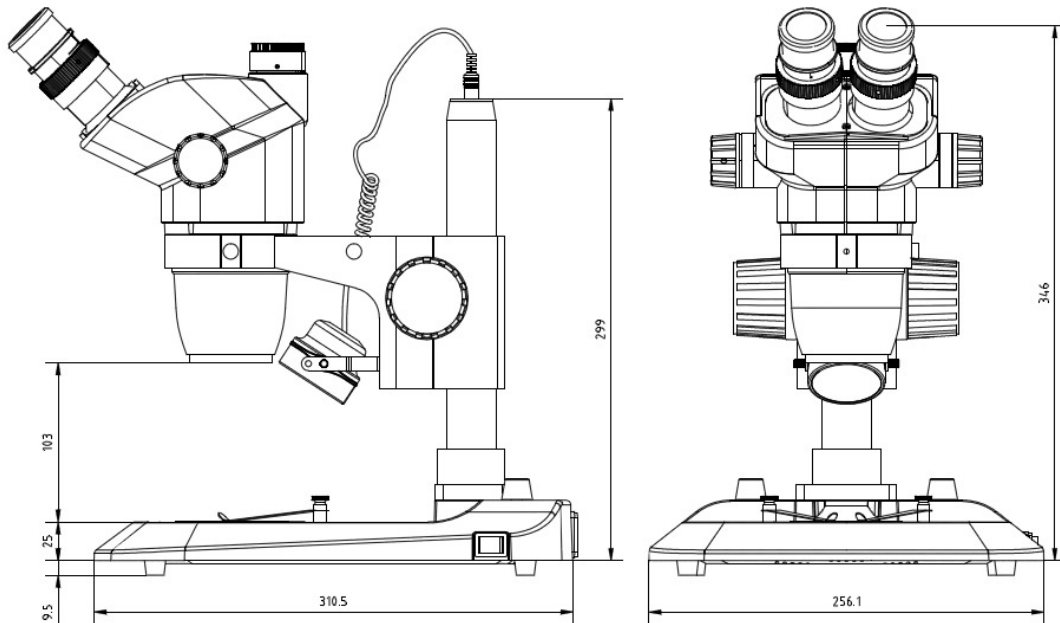Dimensions


NZ.1903-P

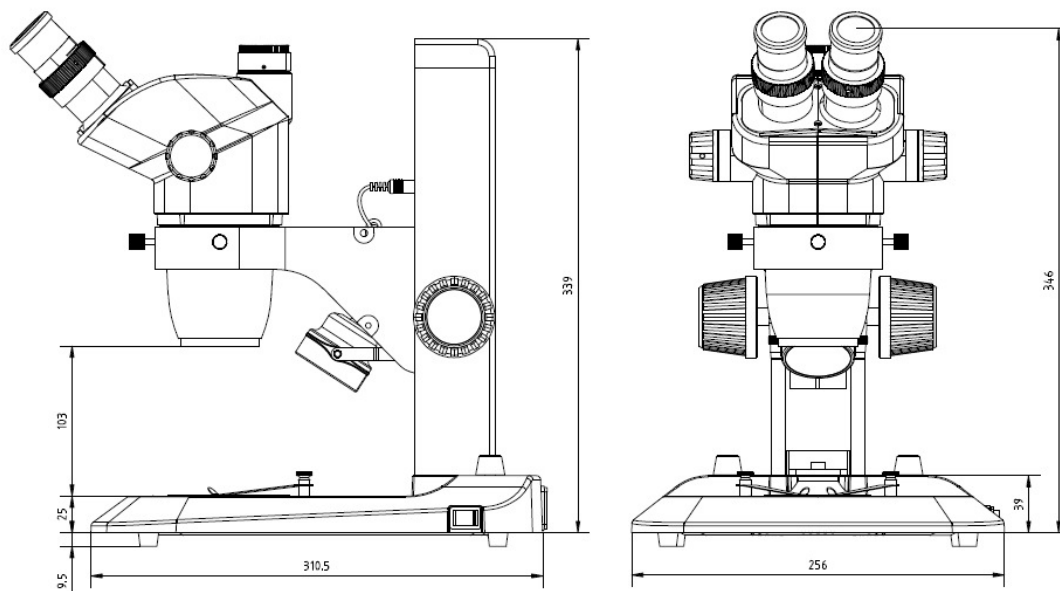

NZ.1903-S

Datasheet v.324112

euromex microscopen

 info@euromex.com
 www.euromex.com

 Papenkamp 20
 6836 BD Arnhem
 Netherlands
 Tel: +31 (0) 26 323 22 11

 Carretera de Barcelona 88, Entresuelo
 Edificio Technomar
 08302 Mataró, Spain
 Tel: +34 717 168 377